

ThetaVAC-2 Vacuum Cold Mounting Machine







地址: 江苏省苏州市吴中区甪直镇龚塘路 228 号 26 栋

Add: Building 26, No. 228, Gongtang Road, Wuzhong District, Suzhou, China Tel: +86-512-65283666 www.trojanchina.com sales@trojanchina.com



Product Introduction:


Trojan's new LCD touch screen programmable vacuum system for cold mounting in epoxy offers excellent pore impregnation in a compact format. This equipment can make epoxy resin penetrate into any porous samples quickly and effectively, resulting in ensured the edge protection of the sample. This equipment is light and compact, external connecting the compressed air source can quickly and effectively pulls a vacuum, and can also be connected to vacuum air pump.

Product Features

- Ability to set the vacuum level, time under vacuum and the number of cycles under vacuum.
- Large Sample Chamber, 10inch (250mm) diameter vacuum chamber
- Easily position samples when pouring under vacuum with electric direction turn table.
- External connecting the compressed air source, easy to operate, safe and reliable (can also be connected to vacuum air pump).
- The inner wall is Teflon coated. Even if the epoxy resin sticks to the inner wall, it is very easy to clean.

地址: 江苏省苏州市吴中区甪直镇龚塘路 228 号 26 栋



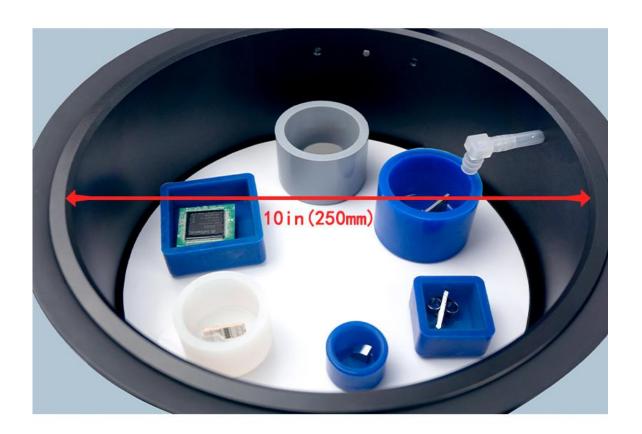

地址: 江苏省苏州市吴中区甪直镇龚塘路 228 号 26 栋


Add: Building 26, No. 228, Gongtang Road, Wuzhong District, Suzhou, China Tel: +86-512-65283666 www.trojanchina.com sales@trojanchina.com



Technical parameter


Electrical	Power Supply	220V, 50Hz, 35W (excluding vacuum pump)		
	Control	5-inch touch screen		
	Timing cycle	Timing range: 0-10min, automatic air extraction /		
		deflation, automatic circulation		
	Program	20 sets		
Vacuum	Vacuum Level	0 ~ - 80kPa, vacuum pump 0 ~ - 90kpa		
	Maximum vacuum flow	75 L/min		
	Air consumption	28-42 L/min		
	Noise	55-60dBA		
Table	Size	The inner diameter is \$\psi\$ 252, and 12 \$\psi\$ 40 samples		
Table		can be placed at the same time		
	Control	Button inching, direction: clockwise / counterclockwise		
Use	Temperature	0-60°C		
environment	Humidity	0-95%RH		
Supply air	Pressure	ϕ 6 gas pipe, the recommended pressure is 0.6 \sim		
		0.7MPa		
External	Voltage	220V~/50Hz		
vacuum pump	Pumping rate	5.4m³/h		
(optional)	Extreme pressure	2pa		
Weight		18.5KG		

Overall dimension

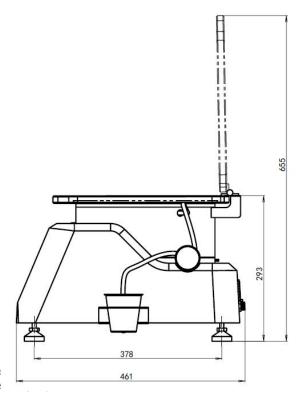

地址: 江苏省苏州市吴中区甪直镇龚塘路 228 号 26 栋 Add: Building 26, No. 228, Gongtang Road, Wuzhong Disti Tel: +86-512-65283666 www.trojanchina.com



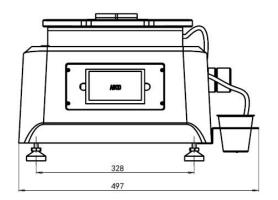

Accessories list

Order No	Name	Specifications	Qty
80.90.100	ThetaVAC-2 Vacuum Cold	86-230V, 50-60Hz	1
	Mounting Machine		
	Enlarged base plate		1
	Seal ring		1
	Disposable paper cup		10 PCS/Pkg
	The stirring rod		10 PCS/Pkg
	Disposable pouring pipe		10 PCS/Pkg
	Epoxy Resin		
	Anti sticking paper		10 PCS/Pkg
	Vacuum pump		1
	Power cord		1
	Instruction manual		1

地址: 江苏省苏州市吴中区甪直镇龚塘路 228 号 26 栋

Add: Building 26, No. 228, Gongtang Road, Wuzhong District, Suzhou, China Tel: +86-512-65283666 www.trojanchina.com sales@trojanchina.com